


CODAN Check Valves



The versatile program
of check valves



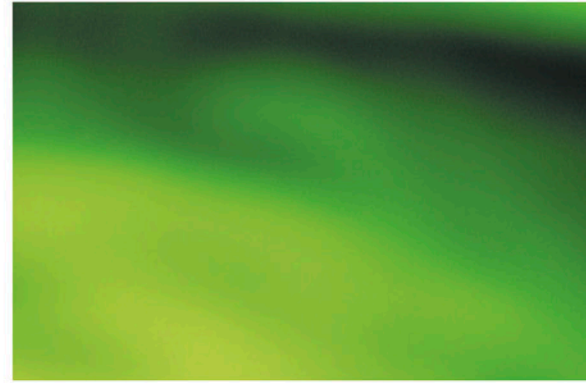

CODAN Check Valves

Integrated safety

CODAN Check Valves are a small detail with maximum importance for patient safety. As a firm or optional integral part in various CODAN products, the check valves reliably prevent reflux or back flow of infusion solution without additional measures and also ensure that blood does not flow back from the patient into the infusion system.

Small, perfect and safe


Top-class engineering in a mini-format

As with all CODAN products, the CODAN Check Valves also manifest decades of experience in research and development.

Precision in production ensures maximum reliability and ease of movement of the CODAN Check Valve. The valve opens upon the slightest conveying pressure. It closes automatically as soon as a counter-pressure results or the conveying pressure has dropped again. Even highly viscous infusions pass through the valves securely and reliably.

The CODAN Check Valves have a particularly low priming volume of 0,15 ml and can be integrated into all infusion systems without any problems.

Be it in standard infusion therapy on wards, in the operating theatre, in intensive care or in anaesthetics – CODAN Check Valves can be used in all areas of everyday hospital work.

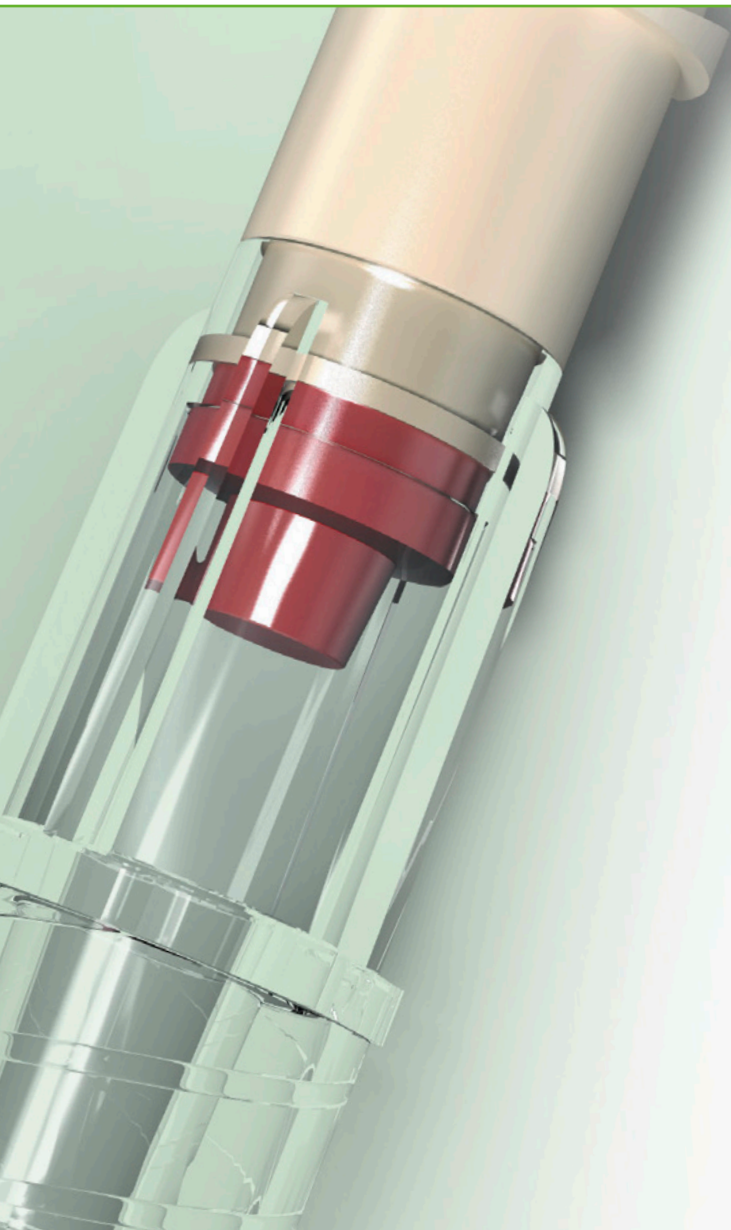

Versatility with precision

For use in many versions


The matching solution can be found for each application. The CODAN Check Valve can be integrated into an existing system as an individual component. In addition, a variety of connection sets, extensions, connectors, completely pre-fabricated systems and special systems with CODAN Check Valves are available. All variants are available with differing tubing diameters. Connectors with CODAN Check Valves are available as double, triple or quadruple connectors.

Connections


CODAN Check Valves in CODAN infusion set connections make handling particularly safe and simple – the less individual components are used, the more quickly and more reliably the infusion can be prepared and administered. Special importance is given to the Y-shaped injection piece in safe and low-risk administration of cytostatics.

Connection variants


The connections provided with a rotating collar which enable secure closure of the Male Luer-Lock without tube torsion are particularly comfortable and user-friendly. Of course CODAN also provides all standard infusion sets with integrated CODAN Check Valves.

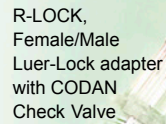


R-LOCK,
Female/Male
Luer-Lock adapter
with CODAN
Check Valve

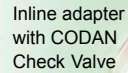


Inline adapter
with CODAN
Check Valve

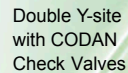


Double Y-site
with CODAN
Check Valves

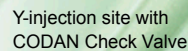


Y-injection site with
CODAN Check Valve

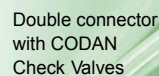

Double connector
with CODAN
Check Valves

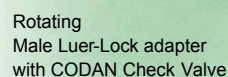

Rotating
Male Luer-Lock adapter
with CODAN Check Valve

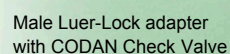

Male Luer-Lock adapter
with CODAN Check Valve



CODAN

CODAN Worldwide



CODAN is known internationally as a manufacturer and supplier of disposable medical transfer systems. The CODAN Companies have more than 1500 employees around the world.

The name CODAN is synonymous with reliability, quality and precision based on the know-how and experience gained from more than 60 years of research and development. Company-owned production facilities and sales companies around the world are a guarantee for efficient production, a tight-knit sales network and a first-class service for customers in the healthcare sector.

CODAN Product range

- Infusion sets
- Transfusion sets
- Extension lines and manifold connectors
- Infusion and transfusion accessories
- Infusion filters and filter systems
- Neonatology/Paediatric products
- Withdrawal, preparation and administration systems
- CODAN CYTO®
- Chemoprotect® products
- Single use syringes
- Invasive blood pressure monitoring systems
- Infusion pumps
- Other CODAN Products

CODAN Companies

CODAN Medizinische Geräte GmbH & Co KG · Deutschland
CODAN pvb Critical Care GmbH · Deutschland
CODAN pvb Medical GmbH · Deutschland
CODAN PORTUGAL, S.A. · Portugal
CODAN 11, S.A. · Portugal
CODAN US Corporation · California · USA
CODAN Inc. · California · USA
CODAN NORGE AS · Norge
CODAN TRIPLUS AB · Sverige
CODAN Limited · Great Britain
CODAN FRANCE Sarl · France
CODAN Medical AG · Schweiz
CODAN ARGUS AG · Schweiz
CODAN BV · Nederland
CODAN s.r.l. · Italia
CODAN Medical GmbH · Österreich
CODAN Medical ApS · Danmark
CODAN DEHA ApS · Danmark

Compliance of the established quality management systems with the provisions of EN ISO 13485 and the Council Directive 93/42/EEC has been certified by the relevant, competent notified bodies:

TÜV SÜD Product Service GmbH

CODAN Medizinische Geräte GmbH & Co KG
23738 Lensahn, Germany
CODAN pvb Critical Care GmbH
85661 Forstinning, Germany
CODAN US Corporation · Santa Ana, CA 92704, USA
CODAN ARGUS AG · 6340 Baar, Switzerland

Presafe Denmark A/S

CODAN Medical ApS · 4970 Rødby, Denmark

Manufacturer

CODAN Medizinische Geräte GmbH & Co KG
Stig Husted-Andersen Straße 11 · 23738 Lensahn, Germany
Tel. +49 (0) 43 63 · 51 11 · Fax +49 (0) 43 63 · 51 12 14
www.codan.de



The decisive connection